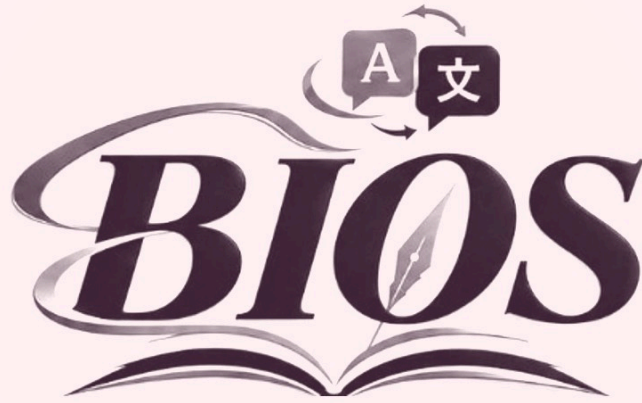
Practical Immunology

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Lab 1 — Safety, Introduction to Serology

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1

Immunity

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EN **Immunity:** body defense against foreign agent.

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The immune response is divided into **cellular** and **humoral** immunity, cellular immunity related to the activity of cells of the immune response, humoral immunity is related to action of some proteins, cytokines and immunoglobulin.

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EN **Lab. Samples:** Blood (Serum – Plasma – RBCs), Urine and cerebrospinal fluid

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EN **Blood:** body fluid in humans and other animals that delivers necessary substances such as nutrients and oxygen to the cells and transports metabolic waste products away from those same cells.

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It is composed of blood cells suspended in blood plasma

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EN Plasma: It is the **yellowish liquid** part of the blood that carries cells and proteins (albumin, globulin and fibrinogen) throughout the body.

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It makes up about **55%** of the body's total blood volume.

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EN Serum: is the fluid and solute component of blood after clotting.

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it is the blood plasma **not including the fibrinogens.**

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Serum includes all proteins not used in blood clotting and all the electrolytes, antibodies, antigens, hormones, and any exogenous substances.

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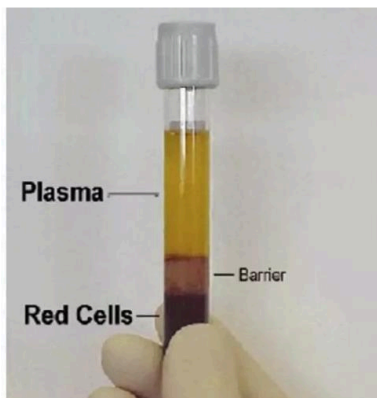
Figure 1. Plasma, barrier and red cells

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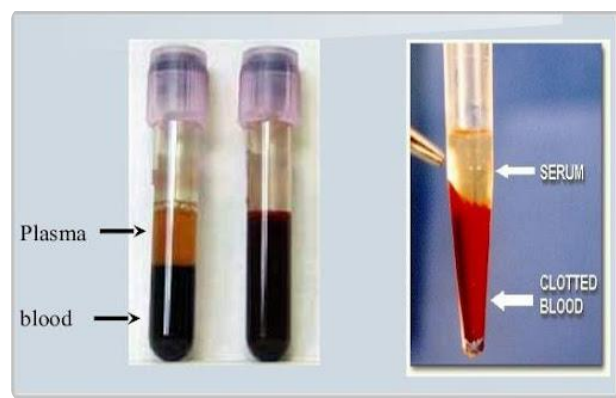
Figure 2. Plasma (left) and serum over clotted blood (right)

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1

EN **Plasma collection** put the blood in tube **with** anticoagulant substance (EDTA, Sodium citrate)

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2

EN **Serum collection** put the blood in tube **without** anticoagulant substance (Plain tube, Gel tube)

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1

EN Collect the blood by vein puncture method.

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2

EN Add the blood into a container with an anticoagulant and then mixed.

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3

EN Add an isotonic solution such as saline to the mixture.

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4

EN Wash the RBCs suspension by centrifugation for **5 min.** at **3000 rpm.**

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5

EN Discard the supernatant and again add saline to the precipitate, mix well and repeat the steps 4 and 5.

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6

EN Re-suspend the precipitate again with saline and standardize it.

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we can preserve the blood in saline for **short time**, and in Alsever's solution for **long time**.

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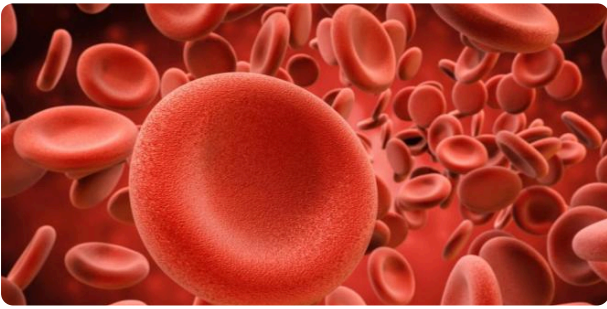
Figure 3. Red blood cells (erythrocytes)

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EN **Saline:** It is normal or isotonic aqueous solutions for sodium chloride in **0.85%** (8.5 gram of NaCl in 1 liter of D.W.).

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- 1 EN Label the tube with the patient's. particulars.

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- 2 EN Put tourniquet on the patient about.
- 3-4'**
- above the venipuncture site.

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- 3 EN Ask patient to form a fist so. veins are more prominent.

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- 4 EN After finding the vein, clean the position

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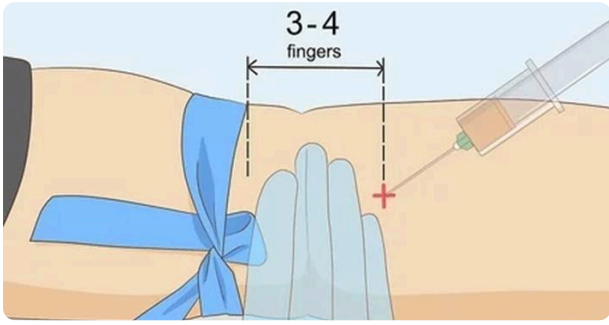
Figure 5. Tourniquet 3-4 fingers above the site

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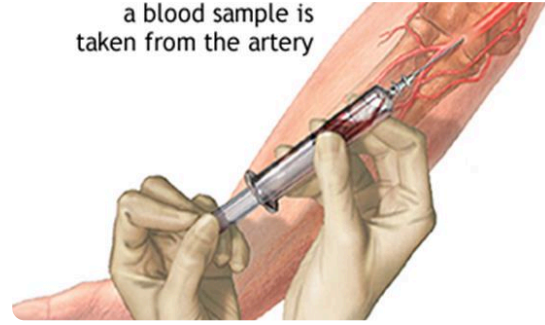
Figure 6. Blood sample taken from the artery after a pulse is found

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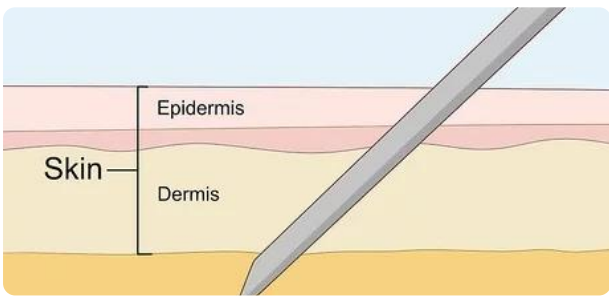
Figure 7. Needle passing through the skin layers

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Safety

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EN

[text cut off in the original file] requires a continuous link between **Three elements:**

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AR

- 1 Source infection: like specimens and contact with them.

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- 2 Method of transmission, and a susceptible host.

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- 3 Procedures used to prevent microorganism transmission which includes:

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EN So, the most safety instructions in immunology laboratory includes:

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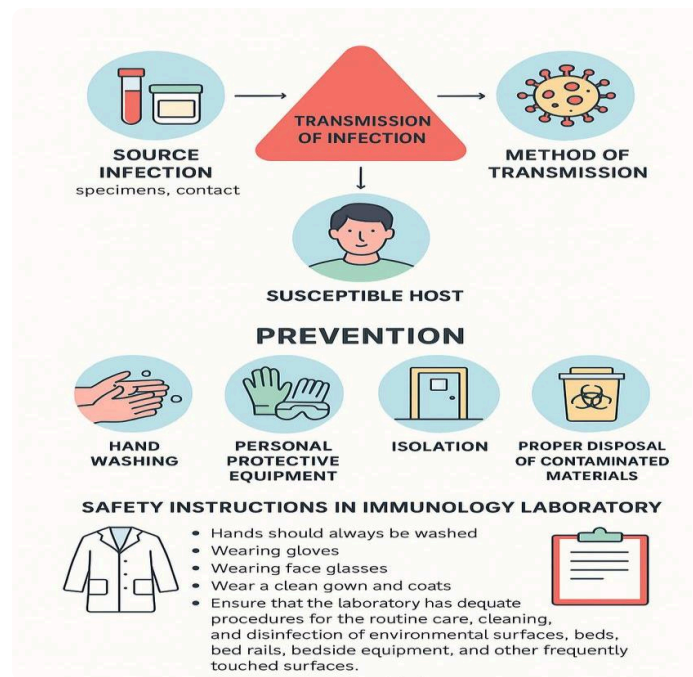
Figure 9. Transmission of infection, prevention and safety instructions

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