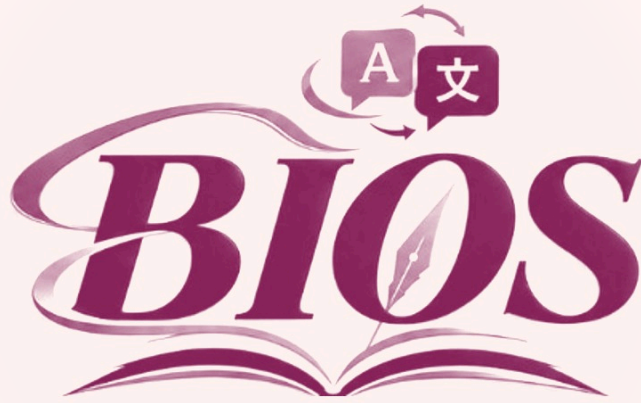
Practical Immunology

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Lab 1 — Safety, Introduction to Serology

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1

Immunity

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EN **Immunity:** body defense against foreign agent. The immune response is divided into **cellular** and **humoral** immunity, cellular immunity related to the activity of cells of the immune response, humoral immunity is related to action of some proteins, cytokines and immunoglobulin.

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EN **Lab. Samples:** Blood (Serum – Plasma – RBCs), Urine and cerebrospinal fluid

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EN **Blood:** body fluid in humans and other animals that delivers necessary substances such as nutrients and oxygen to the cells and transports metabolic waste products away from those same cells. It is composed of blood cells suspended in blood plasma

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EN Plasma: It is the **yellowish liquid** part of the blood that carries cells and proteins (albumin, globulin and fibrinogen) throughout the body. It makes up about **55%** of the body's total blood volume.

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EN Serum: is the fluid and solute component of blood after clotting. it is the blood plasma **not including the fibrinogens**. Serum includes all proteins not used in blood clotting and all the electrolytes, antibodies, antigens, hormones, and any exogenous substances.

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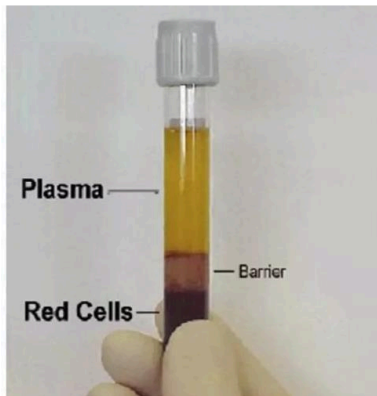
Figure 1. Plasma, barrier and red cells

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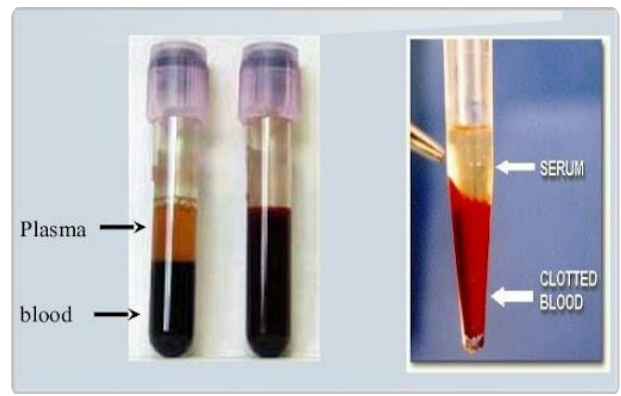
Figure 2. Plasma (left) and serum over clotted blood (right)

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1

EN Plasma collection put the blood in tube **with** anticoagulant substance (EDTA, Sodium citrate)

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2

EN **Serum collection** put the blood in tube **without** anticoagulant substance (Plain tube, Gel tube)

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1

EN Collect the blood by vein puncture method.

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2

EN Add the blood into a container with an anticoagulant and then mixed.

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3

EN Add an isotonic solution such as saline to the mixture.

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4

EN Wash the RBCs suspension by centrifugation for **5 min.** at **3000 rpm.**

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5

EN Discard the supernatant and again add saline to the precipitate, mix well and repeat the steps 4 and 5.

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6

EN Re-suspend the precipitate again with saline and standardize it.

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we can preserve the blood in saline for **short time**, and in Alsever's solution for **long time**.

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Figure 3. Red blood cells (erythrocytes)

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EN

Saline: It is normal or isotonic aqueous solutions for sodium chloride in **0.85%** (8.5 gram of NaCl in 1 liter of D.W.). keep the cell as they are.

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1

EN Label the tube with the patient's particulars.

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2

EN Put tourniquet on the patient about. **3-4'** above the venipuncture site.

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3

EN Ask patient to form a fist so. veins are more prominent.

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4

EN After finding the vein, clean the position

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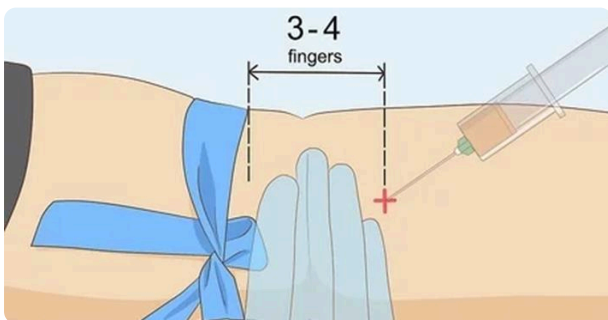
Figure 5. Tourniquet 3-4 fingers above the site

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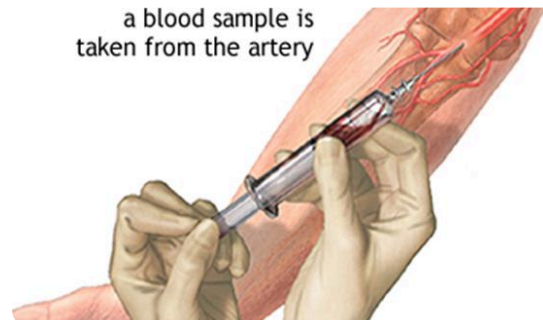
Figure 6. Blood sample taken from the artery after a pulse is found

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Safety

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EN

[text cut off in the original file] requires a continuous link between **Three elements**:

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- 1 Source infection: like specimens and contact with them.

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- 2 Method of transmission, and a susceptible host.

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- 3 Procedures used to prevent microorganism transmission which includes:

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- 3 Isolating highly infective or highly susceptible patients.

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- 4 Properly disposing contaminated materials.

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EN So, the most safety instructions in immunology laboratory includes:

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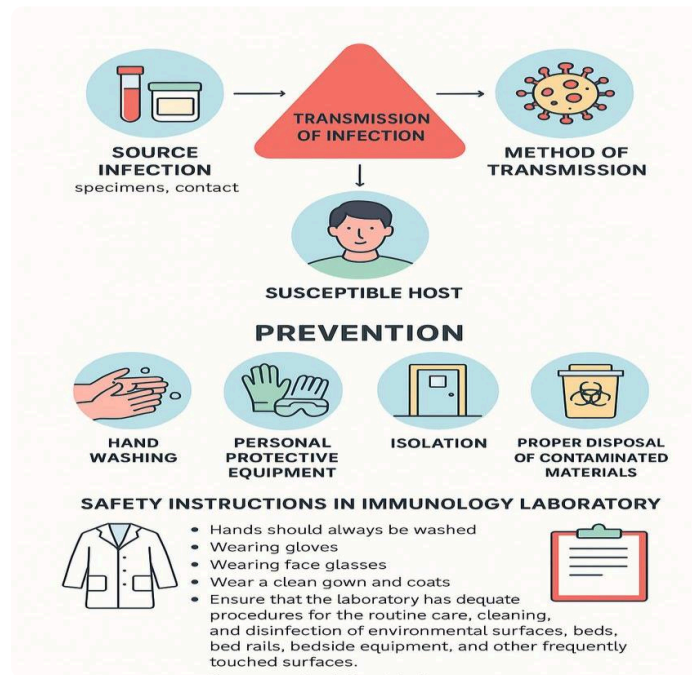
Figure 9. Transmission of infection, prevention and safety instructions

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